

CE  
0123


medicare®  
LIFESENSE®

# User Manual

A1 Blood Pressure Monitor

MD1814 (TMB-1872-A)

Arm Type

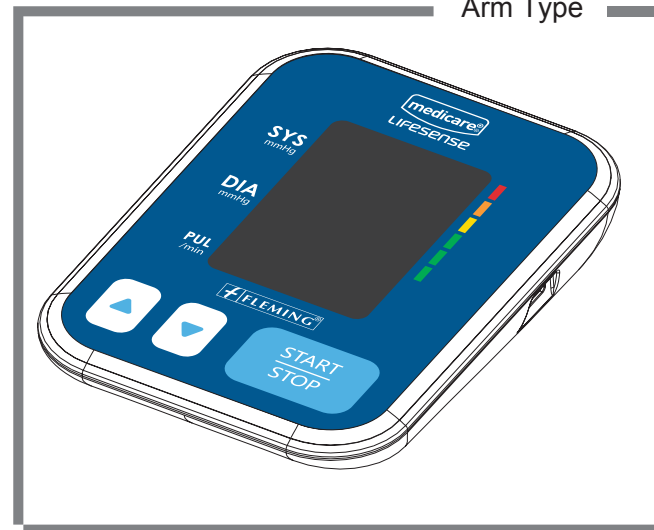

CE  
0123

V1.5.210322

 **GUANGDONG TRANSTEK MEDICAL  
ELECTRONICS CO., LTD.**

Zone A, No. 105, Dongli Road,  
Torch Development District, 528437  
Zhongshan, Guangdong, China

 **MEDICAL DEVICE SAFETY SERVICE GmbH**  
Schiffgraben 41, 30175 Hannover, Germany

 **FLEMING MEDICAL LTD**  
Limerick, V94 K6W4, Ireland  
1800 307777 (IRE) 0800 7314315 (UK)  
[WWW.FLEMINGMEDICAL.IE](http://WWW.FLEMINGMEDICAL.IE)

- Thank you very much for selecting Medicare LifeSense® A1 Blood Pressure Monitor MD1814 (TMB-1872-A).
- Please read the user manual carefully and thoroughly so as to ensure the safe usage of this product. Keep the manual well for further reference in case you have problems.

## Table of Contents

INTRODUCTION .....	2
• General Description	
• Indications for Use	
• Contraindications	
• Measurement Principle	
• Safety Information	
• LCD Display Signal	
• Monitor Components	
• List	
BEFORE YOU START.....	9
• The Choice of Power Supply	
• Installing and Replacing the Batteries	
• Setting Date and Time	
• Select the User	
MEASUREMENT.....	13
• Tie the Cuff	
• Start the Measurement	
DATA MANAGEMENT.....	15
• Recall the Records	
• Delete the Records	
INFORMATION FOR USER .....	18
• Tips for measurement	
• Maintenance	
ABOUT BLOOD PRESSURE .....	20
• What are systolic pressure and diastolic pressure?	
• What is the standard blood pressure classification?	
• Irregular heartbeat detector	
• Why does my blood pressure fluctuate throughout the day?	
• Why do I get a different blood pressure at home compared to the hospital?	
• Is the result the same if measuring on the right arm?	
TROUBLESHOOTING .....	22
SPECIFICATIONS .....	23
AUTHORIZED COMPONENTS.....	24
CONTACT INFORMATION .....	24
COMPLIED STANDARDS LIST.....	25
EMC GUIDANCE.....	26

## ♥ General Description

Thank you for selecting Medicare LifeSense® A1 Blood Pressure Monitor (MD1814 (TMB-1872-A)). The monitor features blood pressure measurement, pulse rate measurement and the result storage. The design provides you with two years of reliable service. Readings taken by the MD1814 (TMB-1872-A) are equivalent to those obtained by a trained observer using the cuff and stethoscope auscultation method. This manual contains important safety and care information, and provides step by step instructions for using the product. Read the manual thoroughly before using the product.

Features:

- 60\*41mm Digital LCD display
- Maximum 60 records per each user
- Measuring during inflation technology

## ♥ Indications for Use

The Medicare LifeSense® Blood Pressure Monitor is a digital monitor intended for use in measuring blood pressure and heartbeat rate with arm circumference ranging from 22cm to 32cm(about 8 $\frac{1}{2}$ "-12 $\frac{1}{2}$ "). It is intended for adult indoor use only.

## ♥ Contraindications










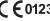


1. The device should not be used by any person who may be suspected of, or is pregnant.
2. The device is not suitable for use on patients with implanted, electrical devices, such as cardiac pacemakers, defibrillators.

## ♥ Measurement Principle

This product uses the Oscillometric Measuring method to detect blood pressure. Before every measurement, the unit establishes a "zero pressure" equivalent to the atmospheric pressure. Then it starts inflating the arm cuff, meanwhile, the unit detects pressure oscillations generated by beat-to-beat pulsatile, which is used to determine the systolic and diastolic pressure, and also pulse rate.

## ♥ Safety Information

The signs below might be in the user manual, labelling or other component. They are the requirement of standard and using.

	Symbol for "THE OPERATION GUIDE MUST BE READ"		Symbol for "TYPE BF APPLIED PARTS"
	Symbol for "MANUFACTURER"		Symbol for "ENVIRONMENT PROTECTION - Electrical waste products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice"
	Symbol for "SERIAL NUMBER"		
	Symbol for "DIRECT CURRENT"		Symbol for "MANUFACTURE DATE"
	Caution: These notes must be observed to prevent any damage to the device.		Symbol for "RECYCLE"
	Symbol for "COMPLIES WITH MDD 93/42/EEC REQUIREMENTS"		Symbol for "MEDICAL DEVICE"
	Symbol for "Authorised Representative in the European Community"		

## ⚠ CAUTION

- \* This device is intended for adult use in homes only.
- \* The device is not suitable for use on neonatal patients, pregnant women, patients with implanted, electronic devices, patients with pre-eclampsia, premature ventricular beats, atrial fibrillation, peripheral, arterial disease and patients undergoing intra-vascular therapy or arterio-venous shunt or people who received a mastectomy. Please consult your doctor prior to using the unit if you suffer from illnesses.
- \* The device is not suitable for measuring the blood pressure of children. Ask your doctor before using it on older children.
- \* The device is not intended for patient transport outside a healthcare facility.
- \* The device is not intended for public use.
- \* This device is intended for non-invasive measuring and monitoring of arterial blood pressure. It is not intended for use on extremities other than the arm or for functions other than obtaining a blood pressure measurement.
- \* Do not confuse self-monitoring with self-diagnosis. This unit allows you to monitor your blood pressure. Do not begin or end medical treatment without asking a physician for treatment advice.
- \* If you are taking medication, consult your physician to determine the most appropriate time to measure your blood pressure. Never change a prescribed medication without consulting your physician.
- \* Do not take any therapeutic measures on the basis of a self measurement. Never alter the dose of a medicine prescribed by a doctor. Consult your doctor if you have any question about your blood pressure.
- \* When the device is used to measure patients who have common arrhythmias such as atrial or ventricular premature beats or atrial fibrillation, the best result may occur with deviation. Please consult your physician about the result.
- \* Don't kink the connection tube during use, otherwise, the cuff pressure may continuously increase which can prevent blood flow and result in harmful injury to the PATIENT.
- \* When using this device, please pay attention to the following situation which may interrupt blood flow and influence blood circulation of the patient, thus cause harmful injury to the patient: connection tubing kinking too frequent and consecutive multiple measurements; the application of the cuff and its pressurization on any arm where intra-vascular access or therapy, or an arterio-venous (A-V) shunt, is present; inflating the cuff on the side of a mastectomy.
- \* Warning: Do not apply the cuff over a wound; otherwise it can cause further injury.
- \* Do not inflate the cuff on the same limb which other monitoring ME equipment is applied around simultaneously, because this could cause temporary loss of function of those simultaneously-used monitoring ME equipment.
- \* On the rare occasion of a fault causing the cuff to remain fully inflated during measurement, open the cuff immediately. Prolonged high pressure (cuff pressure > 300mmHg or constant pressure > 15mmHg for more than 3 minutes) applied to the arm may lead to an ecchymosis.
- \* Please check that operation of the device does not result in prolonged impairment of patient blood circulation.


## ⚠ CAUTION

- \* When measurement, please avoid compression or restriction of the connection tubing.
- \* The device cannot be used with HF surgical equipment at the same time.
- \* The ACCOMPANYING DOCUMENT shall disclose that the SPHYGMOMANOMETER was clinically investigated according to the requirements of ISO 81060-2:2013.
- \* To verify the calibration of the AUTOMATED SPHYGMOMANOMETER, please contact the manufacturer.
- \* This device is contraindicated for any female who may be suspected of, or is pregnant. Besides providing inaccurate readings, the effects of this device on the foetus are unknown.
- \* Too frequent and consecutive measurements could cause disturbances in blood circulation and injuries.
- \* This unit is not suitable for continuous monitoring during medical emergencies or operations. Otherwise, the patient's arm and fingers will become anaesthetic, swollen and even purple due to a lack of blood.
- \* When not in use, store the device with the adapter in a dry room and protect it against extreme moisture, heat, lint, dust and direct sunlight. Never place any heavy objects on the storage case.
- \* This device may be used only for the purpose described in this booklet. The manufacturer cannot be held liable for damage caused by incorrect application.
- \* This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in this booklet.
- \* The equipment is not AP/APG equipment and not suitable for use in the presence of a flammable anaesthetic mixture with air of with oxygen or nitrous oxide.
- \* Warning: No servicing/maintenance while the ME equipment is in use.
- \* The patient is an intended operator.
- \* The patient can measure data and change batteries under normal circumstances and maintain the device and its accessories according to the user manual.
- \* To avoid measurement errors, please avoid the condition of strong electromagnetic field radiated interference signal or electrical fast transient/burst signal.
- \* The blood pressure monitor, its adapter, and the cuff are suitable for use within the patient environment. If you are allergic to polyester, nylon or plastic, please don't use this device.
- \* During use, the patient will be in contact with the cuff. The materials of the cuff have been tested and found to comply with requirements of ISO 10993-5:2009 and ISO 10993-10:2010. It will not cause any potential sensation or irritation reaction.
- \* Adaptor is specified as a part of ME EQUIPMENT.
- \* If you experience discomfort during a measurement, such as pain in the arm or other complaints, press the START/STOP button to release the air immediately from the cuff. Loosen the cuff and remove it from your arm.
- \* If the cuff pressure reaches 40 kPa (300 mmHg), the unit will automatically deflate. Should the cuff not deflate when pressure reaches 40 kPa (300 mmHg), detach the cuff from the arm and press the START/STOP button to stop inflation.
- \* Before use, make sure the device functions safely and is in proper working condition. Check the device, do not use the device if it is damaged in any way. The continuous use of a damaged unit may cause injury, improper results, or serious danger.
- \* Do not wash the cuff in a washing machine or dishwasher!

## ⚠ CAUTION


- \* The service life of the cuff may vary by the frequency of washing, skin condition, and storage state. The typical service life is 10000 times.
- \* It is recommended that the performance should be checked every 2 years and after maintenance and repair, by retesting at least the requirements in limits of the error of the cuff pressure indication and air leakage (testing at least at 50mmHg and 200mmHg).
- \* Please dispose of ACCESSORIES, detachable parts, and the ME EQUIPMENT according to the local guidelines.
- \* Manufacturer will make available on request circuit diagrams, component part lists, descriptions, calibration instructions, etc., to assist to service personnel in parts repair.
- \* The Plug/adaptor plug pins insulate the device from the main supply. Do not position the device in a position where it is difficult to disconnect from the supply mains to safely terminate operation of ME equipment.
- \* The operator shall not touch output of batteries/adaptor and the patient simultaneously.
- \* Cleaning :Dust environment may affect the performance of the unit. Please use the soft cloth to clean the whole unit before and after use. Don't use any abrasive or volatile cleaners.
- \* The device doesn't need to be calibrated within two years of reliable service.
- \* If you have any problems with this device, such as setting up, maintaining or using, please contact the SERVICE PERSONNEL of Fleming Medical Ltd. Don't open or repair the device by yourself in the event of malfunctions. The device must only be serviced, repaired and opened by individuals at authorized sales/service centres.
- \* Please report to Fleming Medical Ltd. if any unexpected operation or events occur.
- \* Keep the unit out of reach of infants, young children or pets to avoid inhalation or swallowing of small parts. It is dangerous or even fatal.
- \* Be careful to strangulation due to cables and hoses, particularly due to excessive length.
- \* At least 30 min required for ME equipment to warm from the minimum storage temperature between uses until it is ready for intended use. At least 30 min required for ME equipment to cool from the maximum storage temperature between uses until it is ready for intended use.
- \* This equipment needs to be installed and put into service in accordance with the information provided in the ACCOMPANYING DOCUMENTS;
- \* Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations, walkie-talkies can affect this equipment and should be kept at least a distance d away from the equipment. The distance d is calculated by the MANUFACTURER from the 80MHz to 5.8 GHz column of Table 4 and Table 9 of IEC 60601-1-2:2014, as appropriate.
- \* Please use ACCESSORIES and detachable parts specified/ authorised by MANUFACTURE. Otherwise, it may cause damage to the unit or danger to the user/patients.
- \* There is no luer lock connectors are used in the construction of tubing, there is a possibility that they might be inadvertently connected to intra-vascular fluid systems, allowing air to be pumped into a blood vessel.
- \* Please use the device under the environment which was provided in the user manual. Otherwise, the performance and lifetime of the device will be impacted and reduced.

## ♥ LCD Display Signal

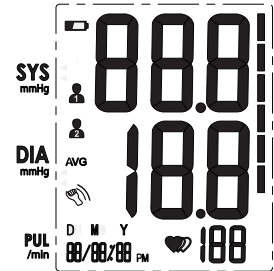


SYMBOL	DESCRIPTION	EXPLANATION
<b>SYS</b>	Systolic blood pressure	High blood pressure
<b>DIA</b>	Diastolic blood pressure	Low blood pressure
<b>PUL /min</b>	Pulse display	Pulse in beats per minute
<b>AVG</b>	Average value	The average value of blood pressure
<b>∩ }</b>	Memory	Indicate it is in the memory mode and which group of memory it is.
<b>Hand</b>	Motion indicator	Motion may result in an inaccurate measurement
<b>mmHg</b>	mmHg	Measurement Unit of the blood pressure (1mmHg=0.133kPa)
<b>⊥ +</b>	Low battery	Batteries are low and need to be replaced
<b>Heart</b>	Irregular heartbeat	Blood pressure monitor is detecting an irregular heartbeat during measurement.
<b> </b>	Blood pressure level indicator	Indicate the blood pressure level
<b>D M Y</b> 00/00/00 PM	Current Time	Year/Month/Day, Hour/Minute
<b>Heart</b>	Heartbeat	Blood pressure monitor is detecting a heartbeat during measurement.
<b>1</b>	User 1	Start measurement for User 1
<b>2</b>	User 2	Start measurement for User 2

## ♥ Monitor Components

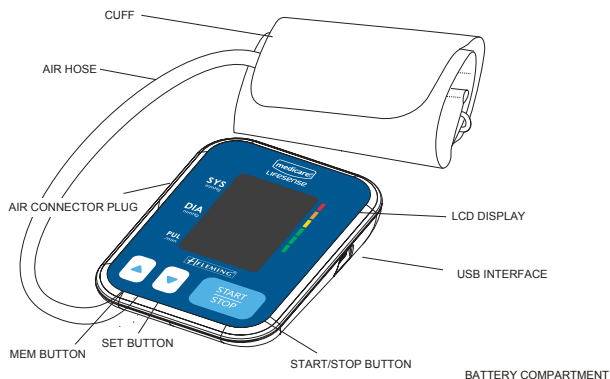

Component list of pressure measuring system

- 1 Cuff
- 2 Air pipe
- 3 PCBA
- 4 Pump
- 5 Valve

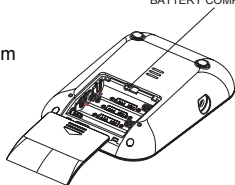

## ♥ List

1. Blood Pressure Monitor MD1814 (TMB-1872-A)
2. Cuff (Type BF applied part) (22~32cm)




3. 4×AAA batteries



4. User manual

## ♥ The Choice of Power Supply

1. Battery powered mode:  
6VDC 4×AAA batteries
2. AC adaptor powered mode: 5V  $\bar{\bar{}}$  1A(Not included)  
Please use the AC adaptor and USB cable just like the following picture:

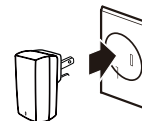


### ⚠ CAUTION


In order to get the best effect and protect your monitor, please use the right battery and special power adaptor which complies with local safety standard.

## ♥ Installing and Replacing the Batteries

- Slide off the battery cover.
- Install the batteries by matching the correct polarity, as shown.
- Replace the battery cover.

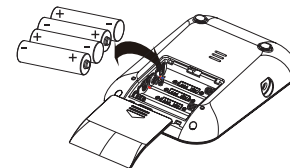

Replace the batteries whenever the below happens

- The  shows
- The display is dim
- The display does not light up




### ⚠ CAUTION

- Do not use new and used batteries together.
- Do not use different types of batteries together.
- Do not dispose the batteries in fire. Batteries may explode or leak.
- Remove batteries if the device is not likely to be used for some time.
- Worn batteries are harmful to the environment. Do not dispose with daily garbage.
- Remove the old batteries from the device following your local recycling guidelines.


## ♥ Setting Date and Time

It is important to set the clock before using your blood pressure monitor, so that a time stamp can be assigned to each record that is stored in the memory. (The setting range of the year :2018—2058 time:24H/12 H)



1. When the monitor is off, press “SET” button shortly it will display the time. Then press “SET” button for 3 seconds to enter the mode for year setting.
2. Press “MEM” button to change the [YEAR]. Each press will increase the numeral by one in a cycling manner.
3. When you get the right year, press “SET” button to set down and turn to next step.

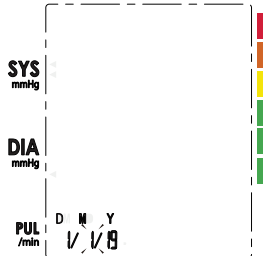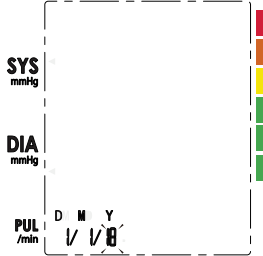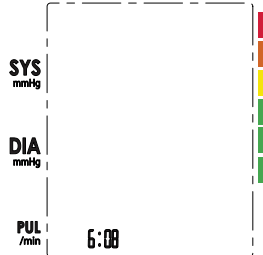



4. Repeat step 2 and 3 to set the [MONTH] and [DAY].

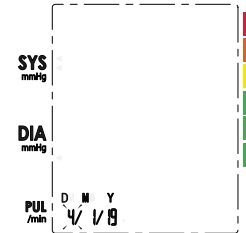

5. Repeat step 2 and 3 to set the [TIME FORMAT] between 12h and 24h.

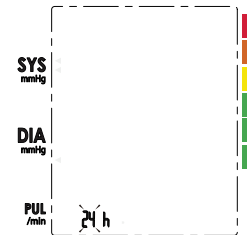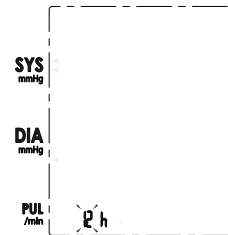


6. Repeat step 2 and 3 to set the [HOUR] and [MINUTE].

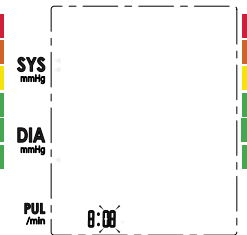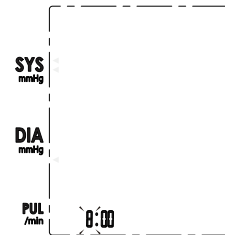


7. After the [MINUTE] is set, the LCD will display “done” first, then display all the settings you have done and then it will turn off.

## ♥ Select the User


- When the monitor is off, press and hold "MEM" button to enter user setting mode. The user ID will blink.

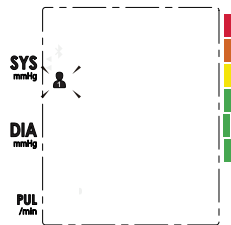

- Then press "MEM" button again, select the user ID between user 1 and user 2.

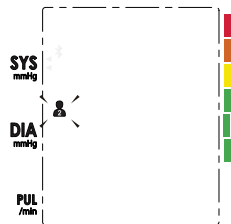


- After selecting the suitable user ID, press "SET" button to confirm. It will display "User ID+ donE" and then turn off.

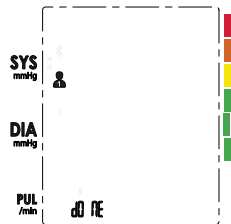


## ♥ Tie the Cuff

- Remove all jewellery, such as watches and bracelets from your left arm. Note: If your doctor has diagnosed you with poor circulation in your left arm, use your right arm.
- Roll or push up your sleeve to expose the skin. Make sure your sleeve is not too tight.
- Hold your arm with your palm facing up and tie the cuff on your upper arm, then position the tube off-centre toward the inner side of arm in line with the little finger. Or position the artery mark  $\Phi$  over the main artery (on the inside of your arm). Note: Locate the main artery by pressing with 2 fingers approximately 2 cm above the bend of your elbow on the inside of your left arm. Identify where the pulse can be felt the strongest. This is your main artery.
- The cuff should be snug but not too tight. You should be able to insert one finger between the cuff and your arm.

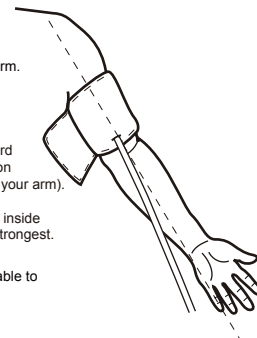

- Sit comfortably with your tested arm resting on a flat surface. Place your elbow on a table so that the cuff is at the same level as your heart. Turn your palm upwards. Sit upright in a chair, and take 5-6 deep breaths.

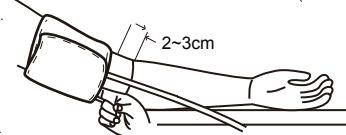

- Helpful tips for Patients, especially for Patients with Hypertension:
  - Rest for 5 minutes before first measurement.
  - Wait at least 3 minutes between measurements. This allows your blood circulation to recover.
  - Take the measurement in a silent room.
  - The patient must relax as much as possible and do not move and talk during the measurement procedure.
  - The cuff should maintain at the same level as the right atrium of the heart.
  - Please sit comfortably. Do not cross your legs and keep your feet flat on the ground.
  - Keep your back against the backrest of the chair.
  - For a meaningful comparison, try to measure under similar conditions. For example, take daily measurements at approximately the same time, on the same arm, or as directed by a physician.






## ♥ Start the Measurement


- When the monitor is off, press “Start/Stop” button to turn on the monitor, and it will finish the whole measurement.




LCD display

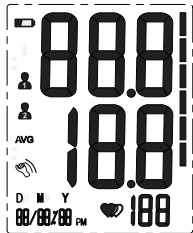


Adjust the zero.

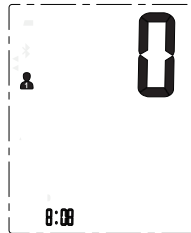

Inflating and measuring.

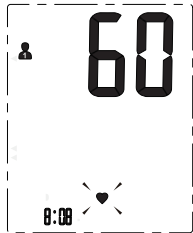

Display and save the results.


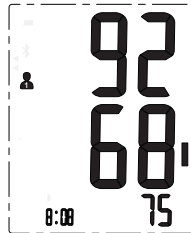


- Press “Start/Stop” button to power off, otherwise it will turn off within 1 minute.



## ♥ Recall the Records

- When the monitor is off, please press “MEM” button to show the recent record. If the records are less than 3 groups, it will display the latest record instead. Take user 1 for example.
- Press “MEM” button or “SET” button to get the record you want.

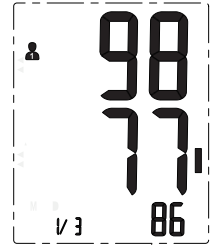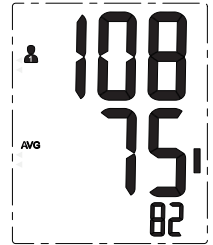

The current No., date and time of the record will be shown alternately.



### ⚠ CAUTION

The most recent record (1) is shown first. Each new measurement is assigned to the first (1) record. All other records are pushed back one digit (e.g., 2 becomes 3, and so on), and the last record (60) is dropped from the list.

## ♥ Delete the Records

If you did not get the correct measurement, you can delete results by following steps below.

### A: Delete one record

1. Hold pressing "set" button for 3 seconds when the monitor is in the memory recall mode, the flash display "dEL one+User ID" will show.
2. Press and hold "Start/stop" button to confirm deleting and the monitor will display "dEL donE+USER ID" and then turn off.



Note: To exit out the delete mode press "Start/Stop" button shortly.





### B: Delete all records

1. Hold pressing "MEM" button for 3 seconds when the monitor is in the memory recall mode, the flash display "dEL ALL+User ID" will show.
2. Press and hold "Start/stop" button to confirm deleting and the monitor will display "dEL donE+USER ID" and then turn off.

Note: To exit out the delete mode press "Start/Stop" button shortly.



3. If there is no record, the right display will show.

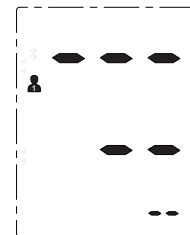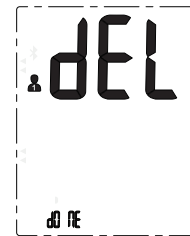

## ♥ Tips for Measurement

Measurements may be inaccurate if taken in the following circumstances.

 <p>Within 1 hour of eating or drinking.</p>	 <p>Immediately after consuming caffeine or smoking.</p>
 <p>Within 20 minutes of taking a bath.</p>	 <p>When talking or moving your fingers.</p>
 <p>When you need to use the toilet.</p>	 <p>In a very cold environment.</p>


## ♥ Maintenance

Please calibrate the blood pressure monitor in specific institute once every two years to ensure the precise measurement. In order to get the best performance, please follow the instructions below for storage.

 <p>Put in a dry place and avoid the sunshine.</p>	 <p>Avoid touching water, clean it with a dry cloth in case.</p>
 <p>Avoid intense shaking and collisions.</p>	 <p>Avoid dusty and unstable temperature environment.</p>
 <p>Using wet cloths to remove dirt.</p>	 <p>Do not attempt to clean the reusable cuff with water and never immerse the cuff in water.</p>


### ♥ What are systolic pressure and diastolic pressure?

When ventricles contract and pump blood out of the heart, the blood pressure reaches its maximum value in the cycle, which is called systolic pressure. When the ventricles relax, the blood pressure reaches its minimum value in the cycle, which is called diastolic pressure.

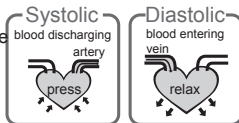

### ♥ What is the standard blood pressure classification?

The blood pressure classification published by World Health Organization (WHO) and International Society of Hypertension (ISH) in 1999 is as follows:

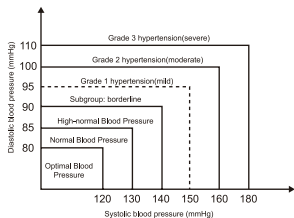

**⚠ CAUTION**  
Only a physician can tell your normal BP range. Please contact a physician if your measuring result falls out of the range. Please note that only a physician can tell whether your blood pressure value has reached a dangerous point.

Level \ Blood Pressure (mm Hg)	Optimal	Normal	High-normal	Mild	Moderate	Severe
SYS	<120	120-129	130-139	140-159	160-179	≥180
DIA	<80	80-84	85-89	90-99	100-109	≥110

### ♥ Irregular Heartbeat Detector

An irregular heartbeat is detected when a heartbeat rhythm varies while the device is measuring systolic pressure and diastolic pressure. During each measurement, blood pressure monitor will keep a record of all the pulse intervals and calculate the average value of them. If there are two or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 25\%$ , or there are four or more pulse intervals, the difference between each interval and the average is more than the average value of  $\pm 15\%$ , then the irregular heartbeat symbol will appear on the display with the measurement result.

**⚠ CAUTION**  
The appearance of the IHB icon indicates that a pulse irregularity consistent with an irregular heart-beat was detected during measurement. Usually this is NOT a cause for concern. However, if the symbol appears often, we recommend you seek medical advice. Please note that the device does not replace a cardiac examination, but serves to detect pulse irregularities at an early stage.

### ♥ Why does my blood pressure fluctuate throughout the day?

1. Individual blood pressure varies multiple times everyday. It is also affected by the way you tie your cuff and your measurement position, so please take the measurement under the same conditions.
2. If the person takes medicine, the pressure will vary more.
3. Wait at least 3 minutes for another measurement.



### ♥ Why do I get a different blood pressure at home compared to the hospital?

The blood pressure is different even throughout the day due to weather, emotion, exercise etc., Also, there is the "white coat" effect, which means blood pressure usually increases in clinical settings.


What you need to pay attention to when you measure your blood pressure at home:  
If the cuff is tied properly.  
If the cuff is too tight or too loose.  
If the cuff is tied on the upper arm.  
If you feel anxious.  
Taking 2-3 deep breaths before beginning will be better for measuring.  
Advice: Relax yourself for 4-5 minutes until you calm down.

### ♥ Is the result the same if measuring on the right arm?

It is ok for both arms, but there will be some different results for different people. We suggest you measure the same arm every time.



This section includes a list of error messages and frequently asked questions for problems you may encounter with your blood pressure monitor. If the product is not operating as you think it should, check here before arranging for servicing.


PROBLEM	SYMPTOM	CHECK THIS	REMEDY
No power	Display will not light up.	Batteries are exhausted.	Replace with new batteries
		Batteries are inserted incorrectly.	Insert the batteries correctly
		AC adaptor is inserted incorrectly.	Insert the AC adaptor tightly
Low batteries	Display is dim or show 	Batteries are low.	Replace with new batteries
Error message	E01 shows	The cuff is too tight or too loose.	Refasten the cuff and then measure again.
	E02 shows	The monitor detected motion, talking, or the pulse is too poor during measurement.	Relax for a moment and then measure again.
	E03 shows	The measurement process does not detect the pulse signal.	Loosen the clothing on the arm and then measure again.
	E04 shows	The treatment of the measurement failed.	Relax for a moment and then measure again.
	EExx, shows on the display.	A calibration error occurred. (XX could be some digital symbol, such as 01, 02, etc., if this similar situation appear, all belongs to calibration error.)	Retake the measurement. If the problem persists, contact the retailer or our customer service department for further assistance. Refer to the warranty for contact information and return instructions.
Warning message	"out " shows	Out of measurement range	Relax for a moment. Refasten the cuff and then measure again. If the problem persists, contact your physician.

<b>Power supply</b>	Battery powered mode: 6VDC 4*AAA batteries AC adaptor powered mode: 5V $\approx$ 1A(not included) (Please only use the recommended AC adaptor model).
<b>Display mode</b>	Digital LCD display V.A.60mm*41mm
<b>Measurement mode</b>	Oscillographic testing mode
<b>Measurement range</b>	Rated cuff pressure: 0mmHg~299mmHg(0kPa ~ 39.9kPa) Measurement pressure: SYS: 60mmHg~230mmHg (8.0kPa~30.7kPa) DIA: 40mmHg~130mmHg (5.3kPa~17.3kPa) Pulse value: (40-199)beat/minute
<b>Accuracy</b>	Pressure: 5°C-40°C within±3mmHg(0.4kPa) Pulse value:±5%
<b>Normal working condition</b>	A temperature range of :+5°C to +40°C A relative humidity range of 15% to 90%, non-condensing, but not requiring a water vapour partial pressure greater than 50 hPa An atmospheric pressure range of : 700 hPa to 1060 hPa
<b>Storage &amp; transportation condition</b>	Temperature:-20°C to +60°C A relative humidity range of ≤ 93%, non-condensing, at a water vapour pressure up to 50hPa
<b>Measurement perimeter of the upper arm</b>	About 22cm~32cm
<b>Net Weight</b>	Approx.162g(Excluding the batteries)
<b>External dimensions</b>	Approx.93mm*130mm*32.5mm
<b>Attachment</b>	user manual,4*AAA batteries
<b>Mode of operation</b>	Continuous operation
<b>Device Classification</b>	Battery Powered Mode: Internally Powered ME Equipment AC Adaptor Powered Mode: Class II ME Equipment
<b>Degree of protection</b>	Type BF applied part
<b>Protection against ingress of water</b>	IP21 It means the device could protected against solid foreign objects of 12.5mm and greater, and protect against vertically falling water drops.
<b>Software Version</b>	A01

WARNING: No modification of this equipment is allowed.

## ♥ Authorized Components

1. Please use the Fleming Medical Ltd. authorised adapter (not included).

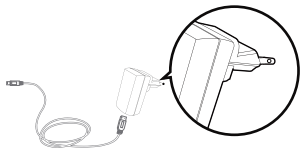

**Adapter**

Type: BLJ06L050100U-V  
 Input:100~240V, 50/60Hz,0.2A MAX  
 Output: 5V = 1000mA

## ♥ Contact Information

For more information about our products, please visit [www.flemingmedical.ie](http://www.flemingmedical.ie). Fleming Medical Ltd. will serve you anytime.

**Distributed By:** Fleming Medical Ltd.  
 Address: Corcanree Business Park, Dock Road Limerick, V94 K6W4, Ireland

Phone (ROI): 1800 307777  
 Phone (UK): 0800 7314315

**Manufactured By:** Guangdong Transtek Medical Electronics Co., Ltd.  
 Address: Zone A, No.105 ,Dongli Road, Torch Development District, 528437 Zhongshan,Guangdong,China

**Authorised EC REP:** Medical Device Safety Service GmbH  
 Address: Schiffgraben 41, 30175 Hannover, Germany

## ♥ Complied Standards List

<b>Risk management</b>	EN ISO 14971:2012 / ISO 14971:2007 Medical devices - Application of risk management to medical devices
<b>Labelling</b>	EN ISO 15223-1:2016 / ISO 15223-1:2016 Medical devices. Symbols to be used with medical device labels, labelling and information to be supplied. Part 1 : General requirements
<b>User manual</b>	EN 1041:2008 +A1:2013 Information supplied by the manufacturer of medical devices
<b>General Requirements for Safety</b>	EN 60601-1:2006+A1:2013+A12:2014 / IEC 60601-1:2005+A1:2012 Medical electrical equipment - Part 1: General requirements for basic safety and essential performance EN 60601-1-11:2015/ IEC 60601-1-11:2015 Medical electrical equipment - Part 1-11: General requirements for basic safety and essential performance - Collateral standard: Requirements for medical electrical equipment and medical electrical systems used in the home healthcare environment
<b>Electromagnetic compatibility</b>	EN 60601-1-2:2015/ IEC 60601-1-2:2014 Medical electrical equipment - Part 1-2: General requirements for basic safety and essential performance - Collateral standard: Electromagnetic disturbances - Requirements and tests
<b>Performance requirements</b>	EN ISO 81060-1:2012 Non-invasive sphygmomanometers - Part 1: Requirements and test methods for non-automated measurement type EN 1060-3:1997+A2:2009 Non-invasive sphygmomanometers - Part 3: Supplementary requirements for electro-mechanical blood pressure measuring systems IEC 80601-2-30:2018 Medical electrical equipment- Part 2-30: Particular requirements for the basic safety and essential performance of automated non-invasive sphygmomanometers
<b>Clinical investigation</b>	EN 1060-4:2004 Non-invasive sphygmomanometers - Part 4: Test procedures to determine the overall system accuracy of automated non-invasive sphygmomanometers ISO 81060-2:2018 Non-invasive sphygmomanometers - Part 2: Clinical validation of intermittent automated measurement type
<b>Usability</b>	EN 60601-1-6:2010+A1:2015/IEC 60601-1-6:2010+A1:2013 Medical electrical equipment - Part 1-6: General requirements for basic safety and essential performance - Collateral standard: Usability IEC 62366-1:2015 Medical devices - Part 1: Application of usability engineering to medical devices
<b>Software life-cycle processes</b>	EN 62304:2006/AC: 2008 / IEC 62304: 2006+A1:2015 Medical device software - Software life-cycle processes
<b>Bio-compatibility</b>	ISO 10993-1:2018 Biological evaluation of medical devices- Part 1: Evaluation and testing within a risk management process ISO 10993-5:2009 Biological evaluation of medical devices - Part 5: Tests for in vitro cytotoxicity ISO 10993-10:2010 Biological evaluation of medical devices - Part 10: Tests for irritation and skin sensitization

## ♥ EMC Guidance

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments

Warning: Don't near active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation."

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment MD1814 (TMB-1872-A), including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

### Technical description:

- All necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the expected service life.
- Guidance and manufacturer's declaration -electromagnetic emissions and Immunity

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class [ B ]
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency
Surge IEC61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV,±2 kV common mode	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV,±2 kV common mode
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % U <sub>r</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % U <sub>r</sub> ; 1 cycle and 70 % U <sub>r</sub> ; 25/30 cycles; Single phase: at 0°.0 % U <sub>r</sub> ; 250/300 cycle	0 % U <sub>r</sub> ; 0.5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % U <sub>r</sub> ; 1 cycle and 70 % U <sub>r</sub> ; 25/30 cycles; Single phase: at 0°. 0 % U <sub>r</sub> ; 250/300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50Hz/60Hz	30 A/m 50Hz/60Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80 % AM at 1 kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80 % AM at 1 kHz
NOTE U <sub>r</sub> is the a.c. mains voltage prior to application of the test level.		

Table 3

Guidance and manufacturer's declaration - electromagnetic Immunity													
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)						
	385	380-390	TETRA 400	Pulse modulation b) 18Hz	1.8	0.3	27						
	450	430-470	GMRS 460, FRS 460	FM c) $\pm$ 5kHz deviation 1kHz sine	2	0.3	28						
	710 745 780	704-787	LTE Band 13, 17	Pulse modulation b) 217Hz	0.2	0.3	9						
	810 870 930							800-960	GSM 800/900, TETRA 800, IDEN 820, CDMA 850, LTE Band 5	Pulse modulation b) 18Hz	2	0.3	28
	1720 1845 1970												
	2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28						
	5240 5500 5785							5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9